



SHELLEY  
MEDICAL  
IMAGING  
TECHNOLOGIES

Sales Office

London, Ontario, Canada

Phone: 1 (519) 690-0874

Fax: 1 (519) 690-0875

Email: sales@simutec.com

Web: www.simutec.com


# Introducer 24F




Introducer 24F lets you introduce your devices (catheters, etc) up to 24F, to our family of models, without leakage.



Connects to 20/26mm silicon tubes.

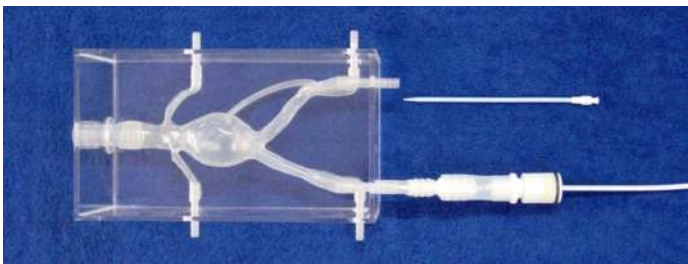

Here the Introducer 24F is in use with a A-S-A-003 model and a 15F catheter going through it.

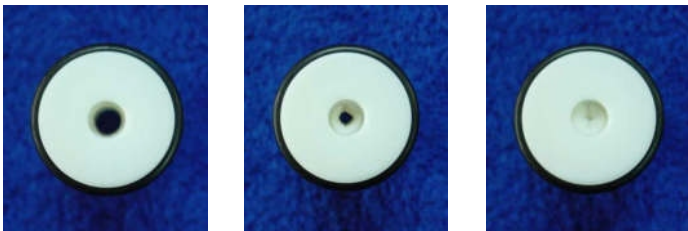

The Introducer 24F can go from wide open to completely closed.



For an even better seal you can cut a glove's finger and fix it on the Introducer (glove not supplied). See also Introducer Luer 10F technical sheet.

ELASTRAT in vitro models are compatible with modern imaging modalities such as digital subtraction angiography, computed tomography and magnetic resonance imaging. Providing the use of an adequate circulating fluid, Doppler techniques can also be performed. The in vitro models transparency to light makes them suitable for video and photographic monitoring.

**ELASTRAT**  
Product Line